

PREMIER elite

Experience the Unmatched Flexibility of
**Operating Space and Deepest Angles on a
Floor Mounted Cath Lab**

CE
2265



Floor Mounted Single Plane Cath Lab


For Cardiovascular, Peripheral Vascular,
Neurovascular Diagnostic and
Interventional Procedures.



L-Arm

Movement

- **±120°** Gantry movement provides unmatched space optimization enabling wide range of Cardiac, Neuro and Peripheral Vascular procedures.
- Enables all Cardiac and Neuro DSA procedures in deeper angles.
- Ease of implantation of all arrhythmia related devices.

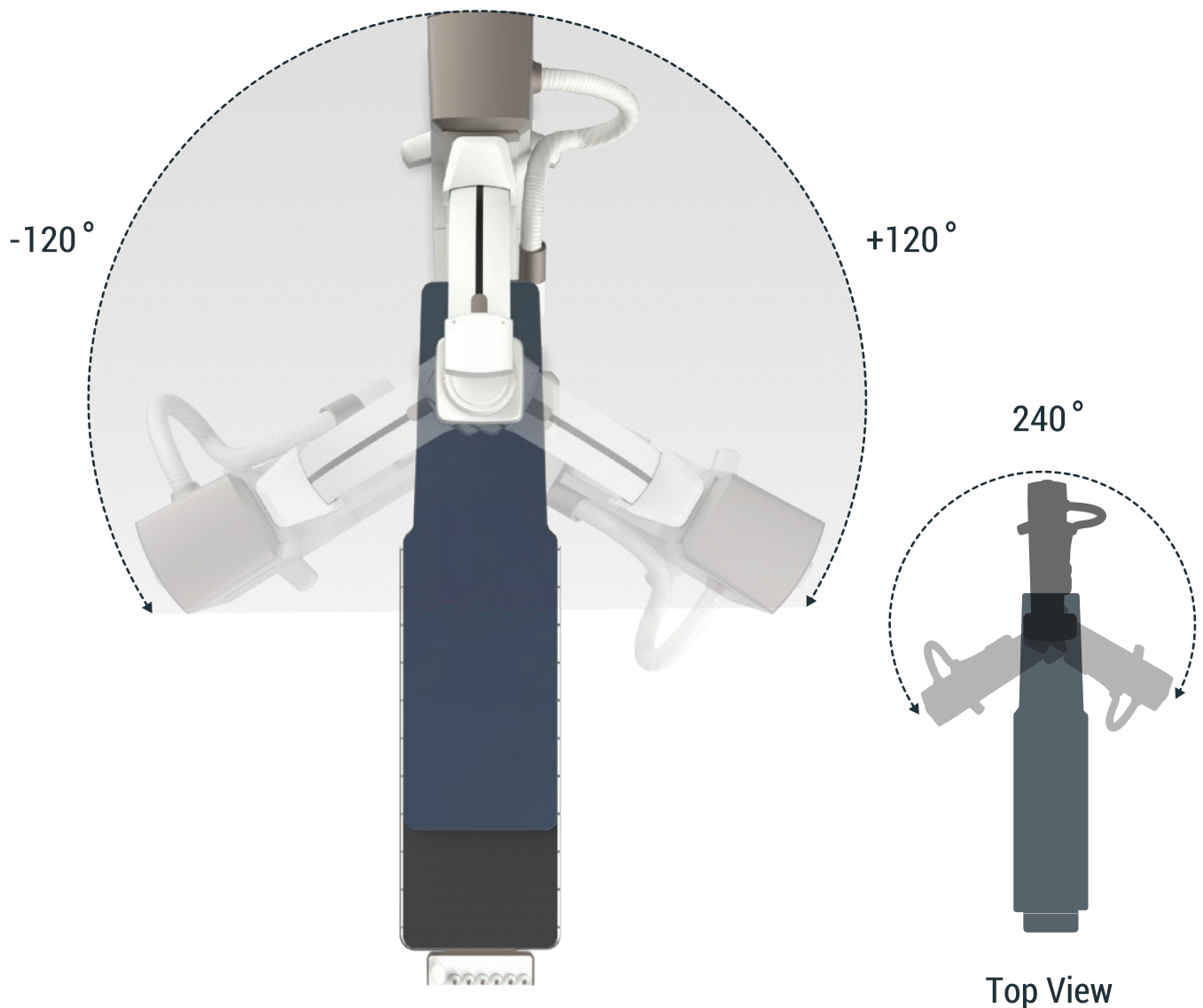

C-Arm Movement


Deeper angles ensure better access to the vasculature and enables better clinical decision making.



Orbital Motion
LAO/RAO
 ± 120 degree



Angular Motion
CRA/CAU
 ± 55 degree

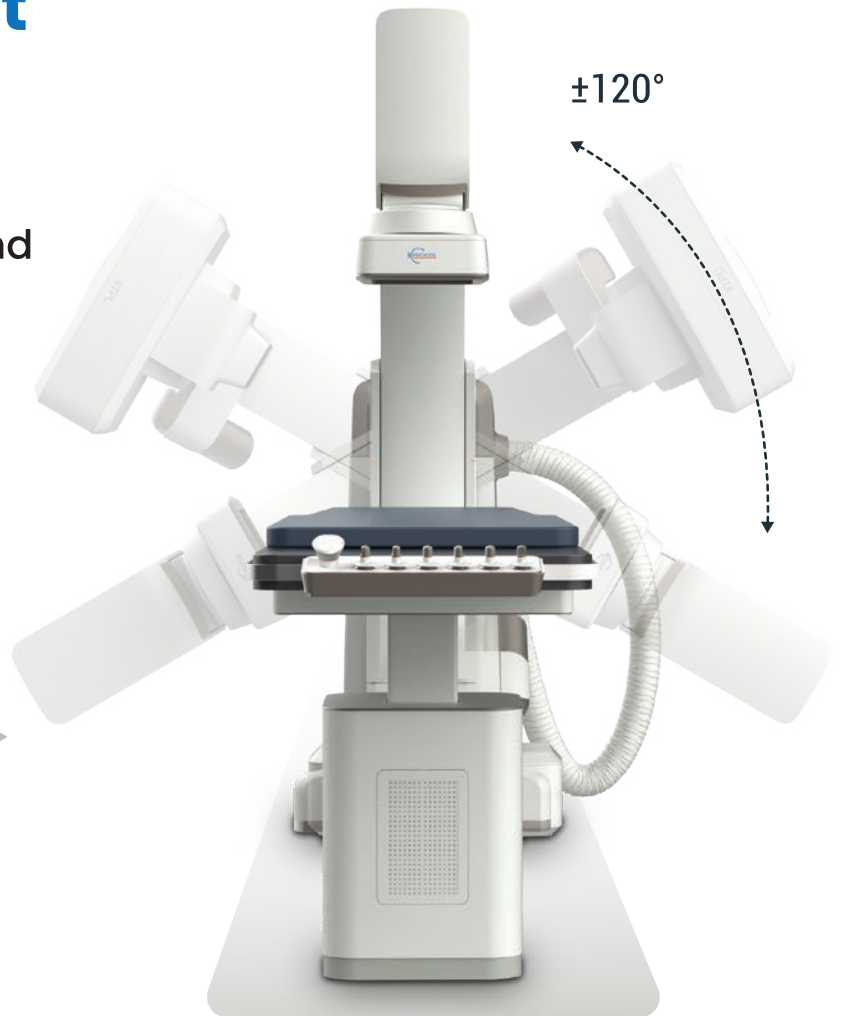

Detector Movement

- Best in class free space for the deepest angles in a floor mounted lab.
- Greater depth in detector movement enables deeper angulations thus helping better anatomical view of Heart.
- Independent Detector Movement as well as synchronized detector movement with L-Arm.

UP/Down 300 mm

Rotation ± 120 degree

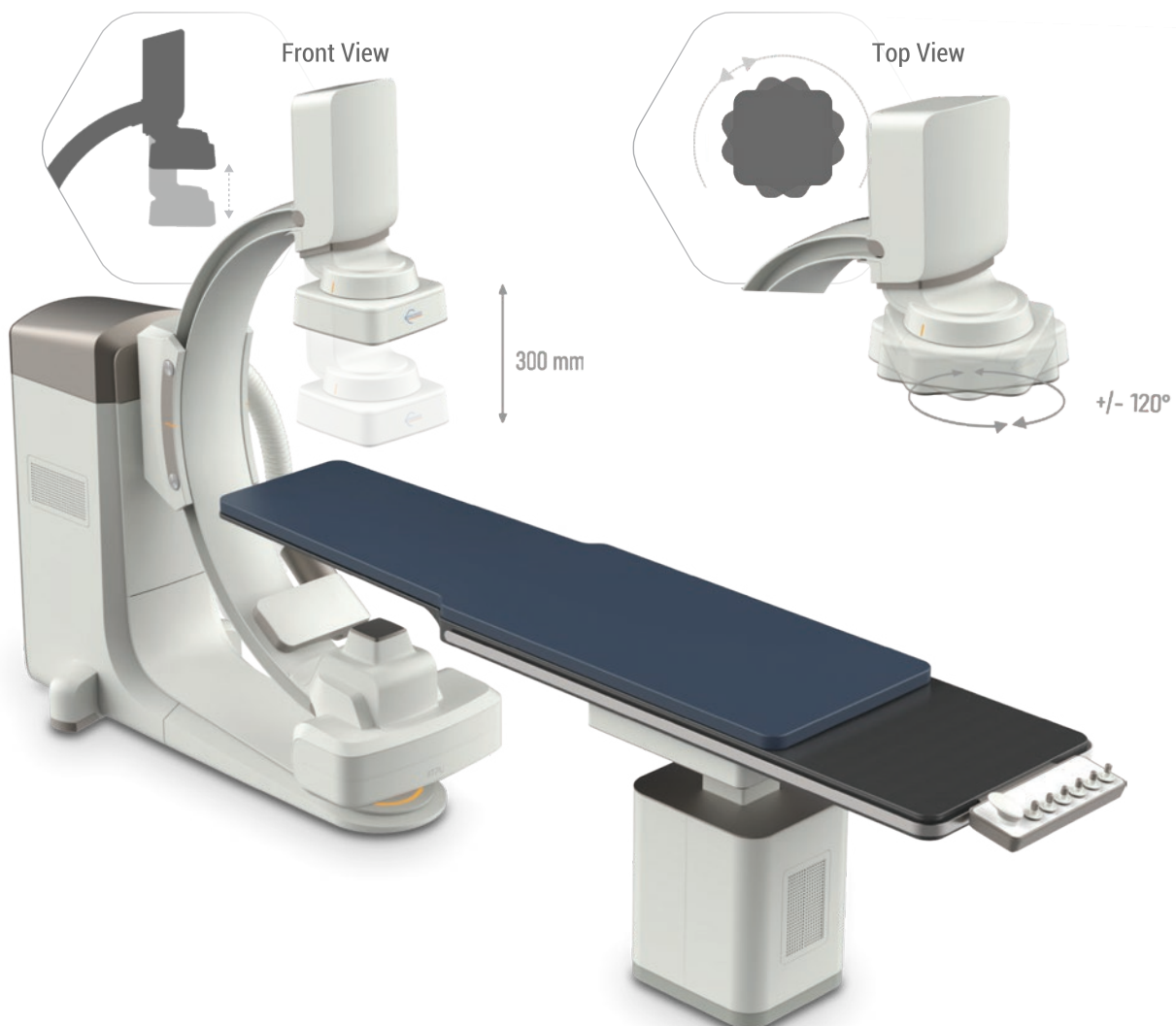


Table Swivel

Table pivot Movement of **$\pm 90^\circ$ Deg**, Provides Greater flexibility to do Vascular Hybrid cases & greater access to the patient for resuscitation during emergencies.

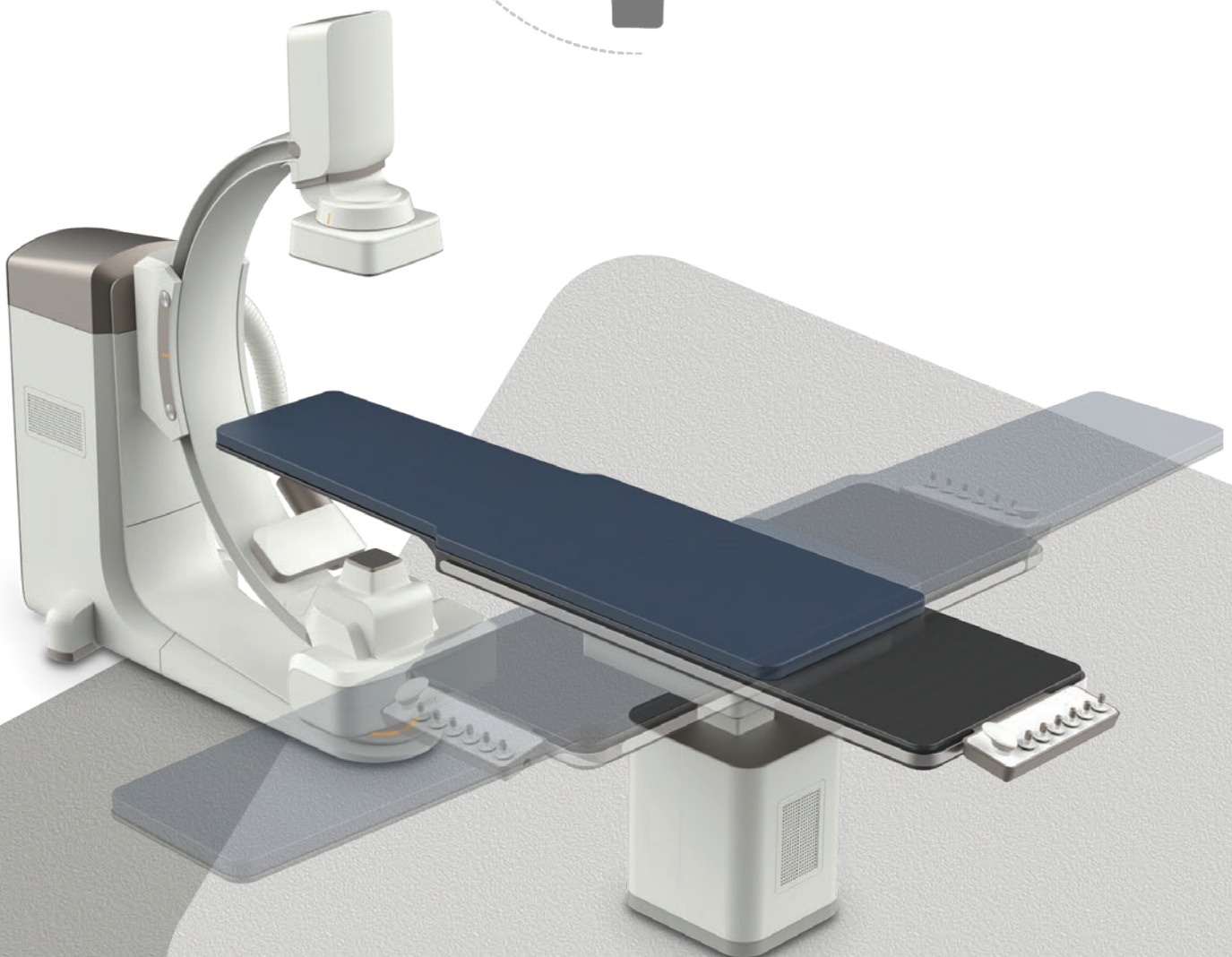
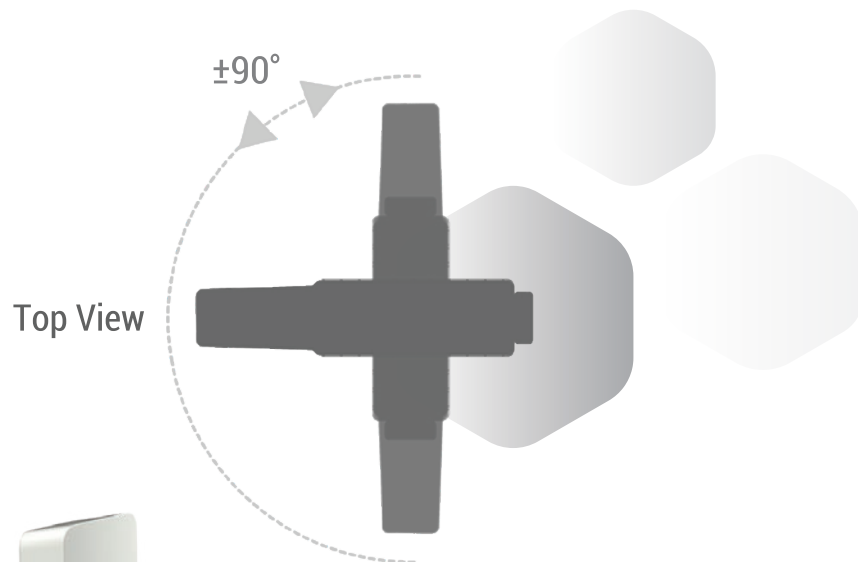


Table Sweep

Head to toe coverage to perform Peripheral Vascular, Neurovascular, Cardiac and lower extremity interventions seamlessly.


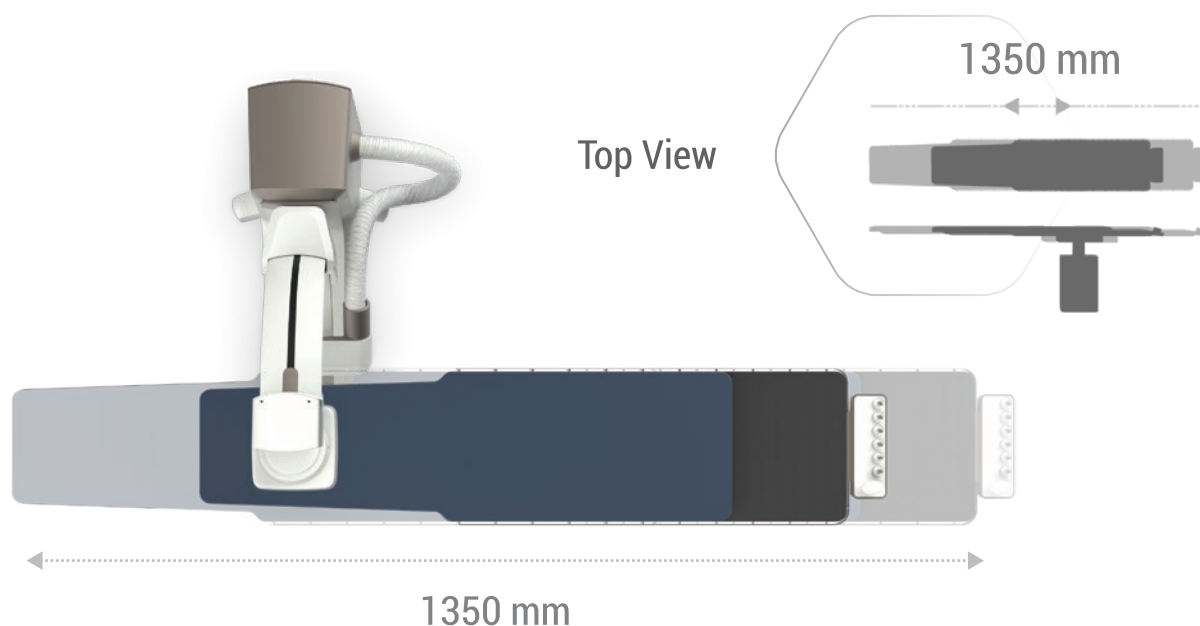



Table Height

Lower table height allows ambulatory patient to sit up & dismount without the need for a patient trolley.

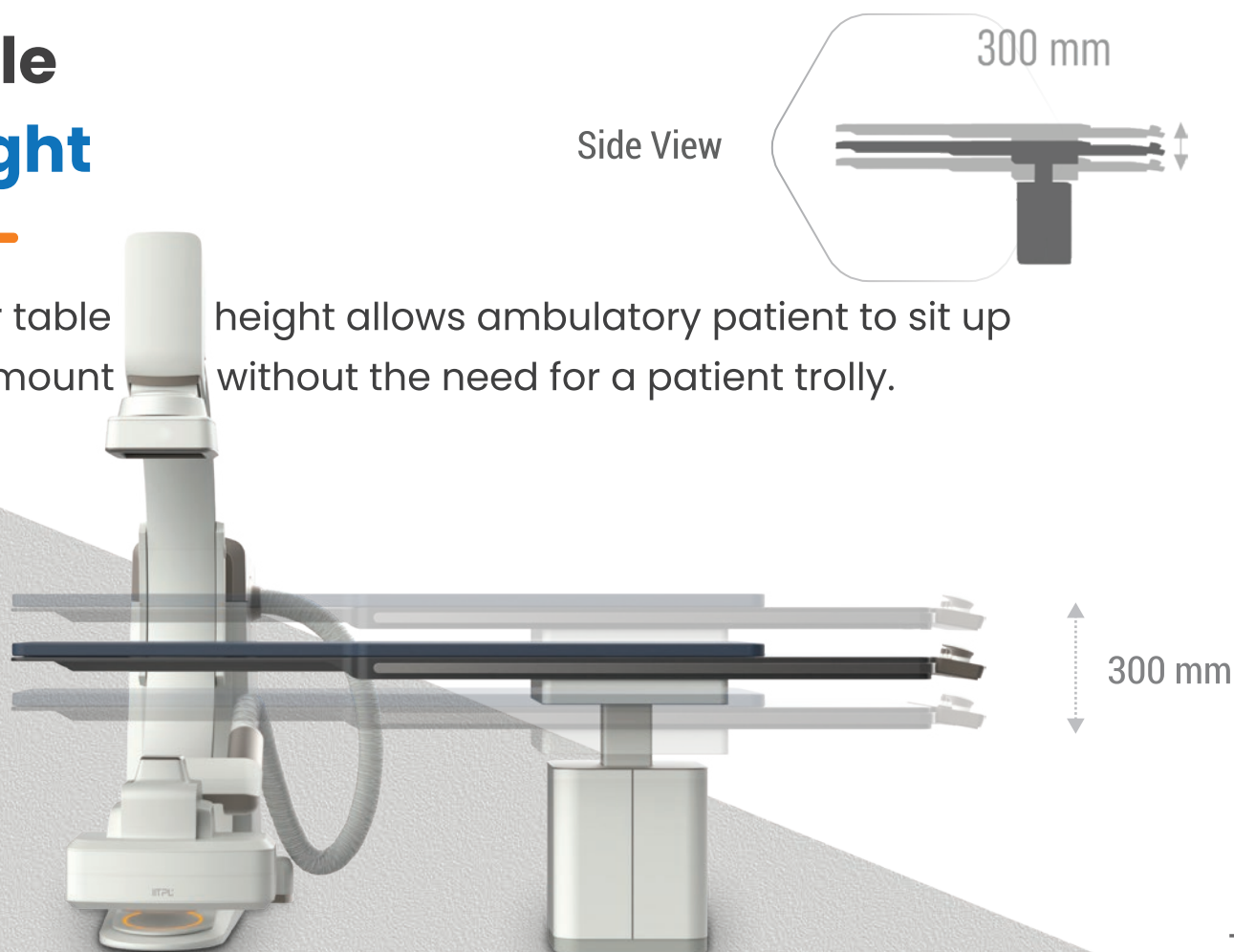

Anticollision Mechanism for Patient Safety

Non-touch based Anticollision
sensors for Detector and L-Arm Side.

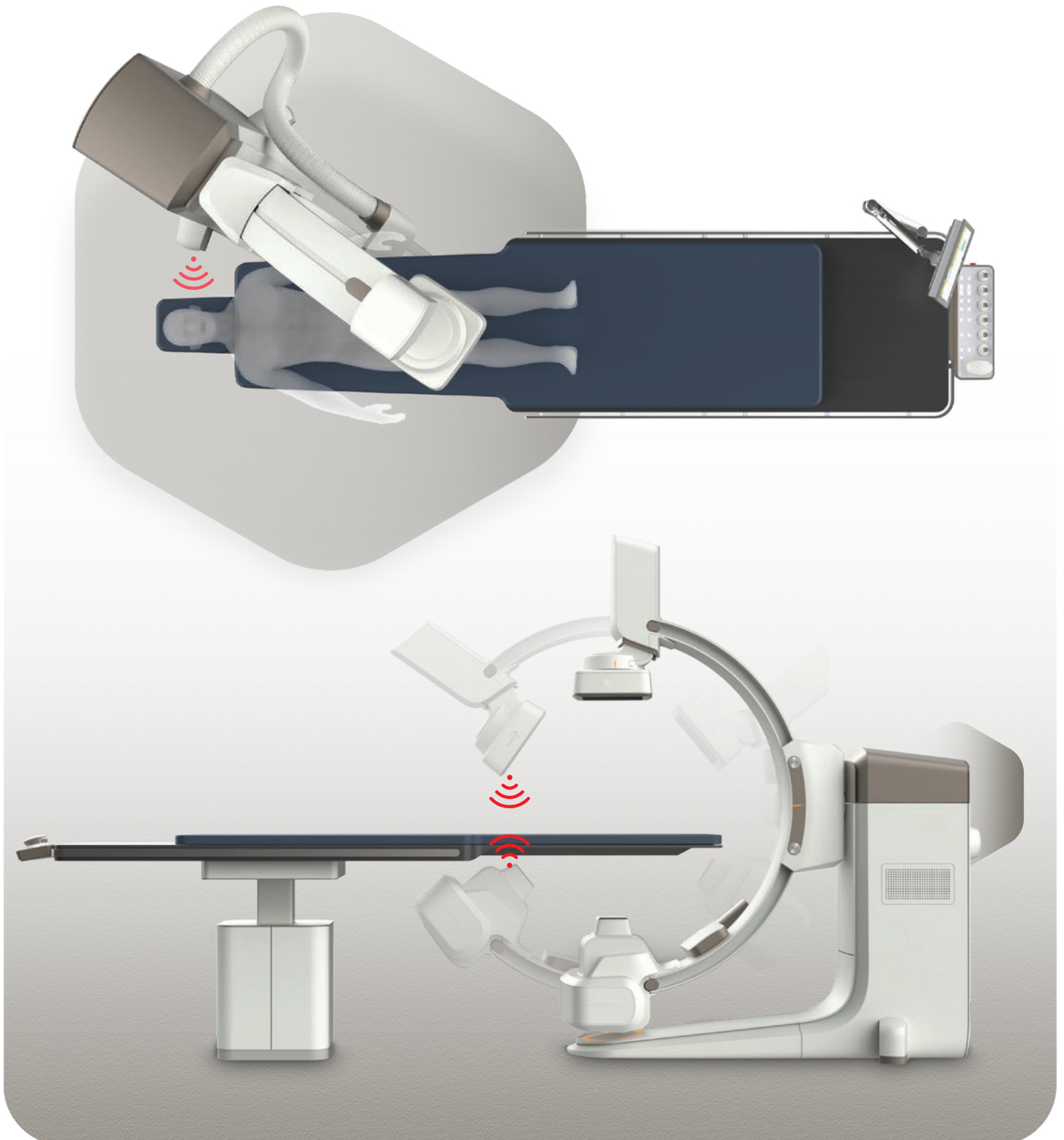

3 MHU Grid Controlled Liquid Metal Bearing Tube

Freedom to carry out lengthy and complex procedures

3 MHU High Capacity Liquid metal bearing tube supports much higher patient throughput and Continuous operation for long hours which improves tube life.



Optimized Dose Management

Grid Controlled Pulsed Fluoroscopy for sharper images and lower radiation.

Reduced Dosage

X-ray tube assembly consists of a metallic case to provide appropriate limits for X-ray.

Enhanced Image Quality

More useful radiation contributes to better image generation.

Liquid metal bearing tube is Virtually silent and delivers near real time imaging.

Smaller Focal Spots (0.4 mm & 0.8 mm)

give much sharper image quality.

100 kW High Frequency X-Ray Generator



- Output power levels of **100 kW**
- **High-Level Radiography**, Continuous Fluoroscopy, Pulsed Fluoroscopy and Cineradiography.
- Better visualization of moving anatomy/OBESE Patients.
- **Higher Switching Frequency** improves image quality & reduces X-Ray dose by increasing X-Ray efficiency.
- **450 KHz** Switching Frequency enables lower high voltage ripples.



Monitor

43" inch

**43" Live Monitor Display
for Sharper & Consistent
Image Quality.**

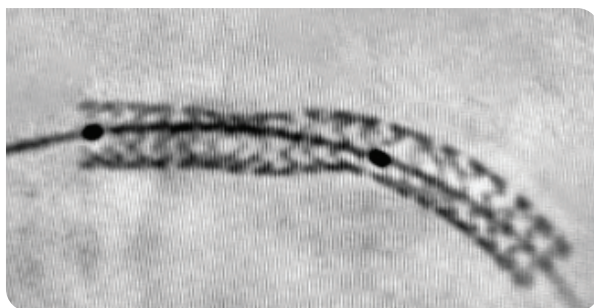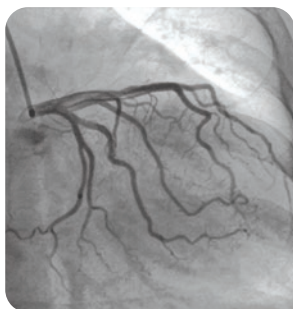

User

Interface

Backlit joysticks, buttons, and keys in the table side module for enhanced user experience.

The table side module includes an **ON/OFF** button for the spot light.



Stent Clarity

Visualise the **entire stent** with **Stent Clarity®** in high resolution.

Optional Features

AINSTA QCA

World's first AI – enabled vessel segmentation analysis, delivering results in just **2 seconds**.




Virtual Fractional Flow Reserve

vFFR is a less invasive angio- derived FFR calculation method. Without the need of adenosine and pressure wire achieving **94% correlation** with invasive FFR measurements to enhance procedural efficiency and patient care.




QCA

Quantitative Coronary Analysis accurately determines vessel diameter and quantify stenosis, which is crucial for selecting the optimal stent dimensions to suit the original vessel diameter and lesion length.

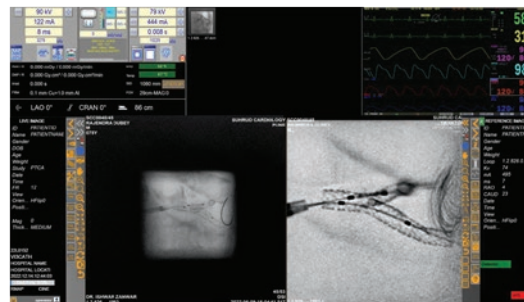

AI-MAG

AI guided advanced magnification and superior resolution with **zero additional radiation exposure** and **zero additional contrast**.

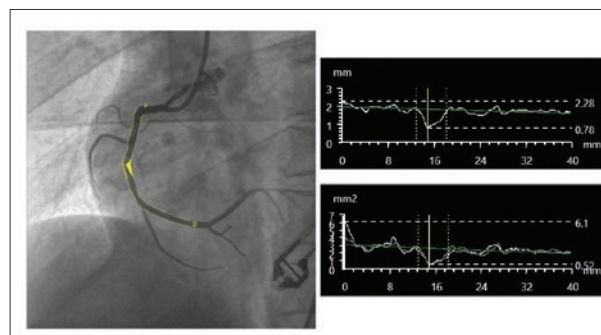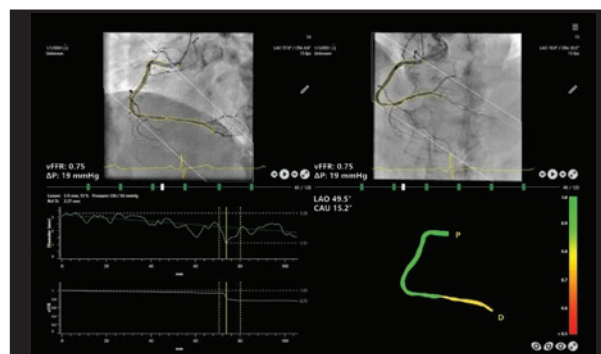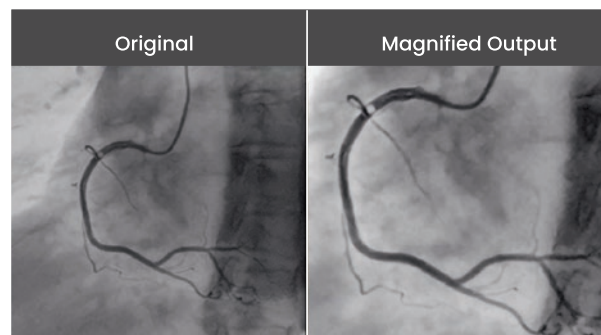


Optional Features


IV-LINQ (Co-registration of IVUS/OCT data with angiography)

Co-register real time between **angiography** and **IVUS/OCT** with **IV-LINQ** to obtain a detailed view of the lesion and its exact location in the coronary tree.

IVUS Co-registration

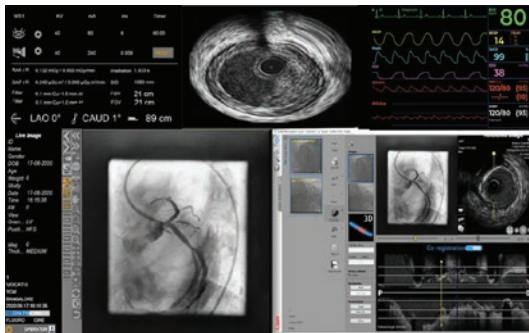

OCT Co-registration

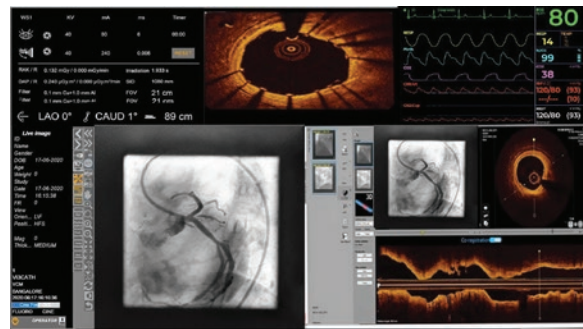

Advanced Ventricular Assessment

Enhances ventricular assessment with **LVA/RVA Analysis** for **automated workflows** to precisely calculate volumes, analyze wall motion and assess ejection percentage.

Precise Peripheral Vessel Analysis

From the **abdominal aorta to carotid, renal, iliac and femoral arteries**, the workflow offers unparalleled precision to improve procedural efficiency and accuracy.

Integrated Ultrasound

With combined precision and **ultrasound guidance, including vascular Doppler studies**, perform guided arterial punctures and transthoracic echocardiography during procedures.

Elevated TAVI Planning

(TAVI) planning before the procedure with the **3mensio Structural Heart software** that enables in elevating procedural success rates with detailed visualization and meticulous analysis.

Specification



Gantry

C Arm Speed	: Up to 25 degree/Sec
ISO Center Height	: 1110 mm
Min. SID	: 960 mm
Max. SID	: 1260 mm
SID Travel	: 300 mm
Max Free Space	: 930 mm
L-Movement	: ±120 degree, 3 side of table
Home Position L-arm	: Head End
Speed	: Up to 15 degree/Sec
2 Way Anti-collision	: Detector side and L-Arm side



Table

Table top	: Carbon Fiber
Mattress	: Memory foam
Max. patient weight	: 200 Kg+50 kg CPR
Min. height	: 780 mm
Vertical stroke	: 300 mm
Longitudinal stroke	: 1350 mm
Table lateral stroke	: 300 mm
Table pivot Movement	: ±90 degree



Detector:

Flat Panel	: 12" diagonal
Pixel Size	: 154 µm
Image Matrix	: 1344x1344
Dynamic range	: 16-bit
DQE	: 77%



Tube:

3 MHU Grid Controlled
Liquid Metal Bearing tube



Generator:

100 kW (High frequency
X-Ray generator) 40-125 kV



Display Monitor:

43"or 55" Live Monitor Display.

43" 4K Industrial grade reference monitors in console room.

About us

Innvolution Healthcare is a group focused on interventional cardiology with a comprehensive portfolio of drug-eluting stents, balloon catheters, and the award-winning and fastest-growing Cath Lab Pinnacle & Premier. Innvolution Group comprises **two** companies



Innvolution Healthcare Private Limited ("IHPL"), which manufactures and distributes cardiac care consumables.



Innovation Imaging Technologies Private Limited ("IITPL"), which manufactures Cath Labs.

The company has three manufacturing facilities in India located in **Bengaluru, Vizag, and Jaipur.**

Innvolution is the fastest-growing and Asia's most awarded Cath Lab manufacturing company. The promoters of the company have a vision of creating truly differentiated solutions in Medical Imaging and offering solutions that make cardiac care more effective and accessible.

16* Awards

16 National & International Awards. Asia's Most Awarded Cath Lab Company.



93%*

93% customer satisfaction across 400+ installations in India.



48* Hrs.

Less than 48 Hrs of downtime at any installation since inception.



* Data as of Oct 2024



reddot design
award 2018



GOOD DESIGN
AWARD 2019



National Innovation Award
by President of India,
TDB.



World's Best
Emerging Brand 2023





CE
2265



Address:

121F, Bommasandra Industrial Area, Phase 1, Hulimangla Road,
Hosur Main Road, Electronic City, Bangalore 560 099, Karnataka, India

Phone:

+91 6364 882 807

Mail:

info@innvolution.com

Website:

www.innvolution.com



SAL-ST-010, Rev 02